FOOTBEAT[®]

INSTRUCTIONS FOR USE

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All rights reserved. All trademarks and registered trademarks are the property of their respective owners.

This product is covered by one or more of the following U.S. Patents: 7,909,783, 8,246,556, 9,283,139, and 9,439,828 and various pending U.S. patent applications. Purchase of this package does not act to grant a use license under these patents in any jurisdiction outside of the United States. Additional patents pending.

Information and views expressed in this document, including URL and other internet web site references, may change without notice.

Manufactured by:
AVEX LLC
Grand Junction, CO

Footbeat Support

Toll Free: 1-844-388-AVEX

Hours: 8:00 A.M. – 5:00 P.M. MST

Email: support@footbeat.com

About This Manual

This manual describes the Footbeat circulation enhancement and recovery system ("Footbeat") and provides instructions for its use. Please read the entire manual carefully before use. This manual contains information, contraindications and warnings that must be observed to ensure safe operation and maintain Footbeat in a safe condition.

The Footbeat System

START RECOVERING. Footbeat is the natural, non-invasive solution for enhancing circulation throughout the body when you are seated, immobile or lying down. Footbeat's patented micro-mobile compression technology, speeds recovery and improves quality of life by compressing the veins of the foot and naturally pumping blood up the legs. This unique action duplicates the circulatory benefits received from being active or walking.


Technical Description

Footbeat provides therapeutic massage through a cyclic compression to the arch of the foot when the user is in a non-weight bearing position. The Footbeat insole contains a pad that presses into the arch of the foot approximately every 35 seconds. The pad holds for 2 seconds and then fully retracts. When a user walks or runs, the Footbeat system senses this activity and pauses. Upon return to a seated or recumbent position, the Footbeat system will reactivate the compression cycle.

What's In The Box

1. Quick Start Guide
2. Insoles with engines
3. Socks (2 pairs)
4. Tote bag
5. Moccs
6. Instructions For Use (IFU)
7. Remote
8. Charging cable

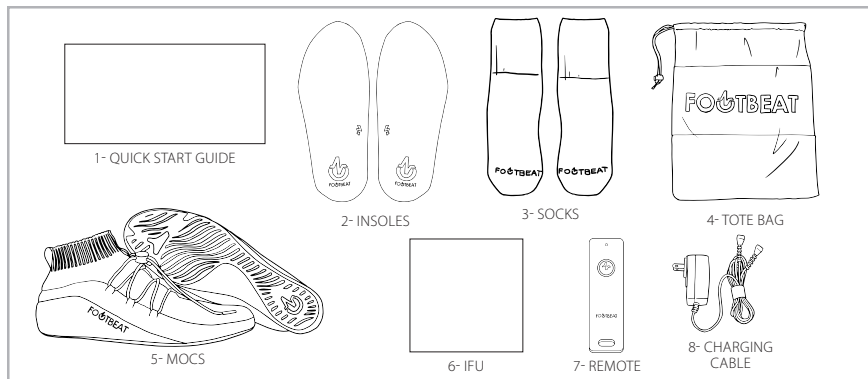Intended Use

Footbeat provides a therapeutic foot massage that provides relief from aches and pains.

NOTE – The instructions are for information only and should not be used for diagnosis, prevention or treatment of medical conditions. Footbeat is not a prescriptive therapy for any disease or illness. Footbeat is non-invasive and not pharmacological. There are no known side effects for Footbeat. Please consult your healthcare provider regarding the duration and frequency that Footbeat should be used to achieve the maximum benefits from Footbeat.

IMPORTANT

Moccs should not be worn outdoors.

Contraindications

Unless directed by a physician, DO NOT use Footbeat if you have:

- a. Recently had a blood clot or are under the care of a physician for Deep Vein Thrombosis (DVT)
- b. Pulmonary edema
- c. Pulmonary embolism (PE)
- d. Foot ulceration, wounds or lesion at or near the plantar venous plexus on the bottom of the foot

Warnings

1. DO NOT use in conjunction with compression stockings or compression sleeves worn around the feet, ankles or legs.
2. DO NOT wear or use while operating any type of vehicle.
3. DO NOT use with ice or cold therapy devices.
4. Failure to observe the following procedures may result in battery leakage, heat, bursting, and serious personal injury.
 - a. Do not remove any batteries from product.
 - b. Never dispose of batteries in fire.
 - c. Keep batteries away from combustible materials.
 - d. Do not incinerate or heat batteries.
 - e. Do not disassemble, alter, or reconstruct product.
 - f. Batteries contain a strong colorless electrolyte, which may cause permanent loss of eyesight and injury when in direct contact with eyes or skin. Immediately flush with clean water when in contact with leaked electrolyte and consult a physician.

To insure proper safety and optimum Footbeat operation, observe the following warnings:

1. There are no serviceable parts. To avoid the risk of electrical shock, do not attempt to open or modify the insoles, charger or remote.
2. To avoid electrical shock, product malfunction or damage, never operate with mechanical damage, or if the unit is otherwise not fully operational.
3. To avoid the risk of electrical shock, unplug the device when in use.
4. Never operate with any AC adapter other than that provided by AVEX LLC. Use of other adapters may result in electrical shock and will void the warranty.

Precautions

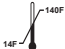
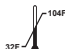
1. Anyone under the age of 18 years old should consult a healthcare provider before using Footbeat.
2. Use the device on the foot only and per instructions.
3. Not intended for use during athletic activities.
4. Because Footbeat is in sleep mode when walking or running and won't compress the foot, it should not be worn for extensive walking. Do not run or jump while wearing Footbeat.
5. Other conditions are presented in this IFU in appropriate locations, identified with the caution symbol.

IMPORTANT

If you have any questions about any existing health conditions, consult your Health Care Provider prior to using Footbeat.

Symbols and Explanations

	Electrical Caution	Refer to the warnings to avoid electrical shock.
	Caution	A hazardous situation, which if not avoided, could result in injury.
	Consult Instructions For Use	Read Instructions For Use for care, storage, operation and warnings.
	Keep Dry	Follow all recommended washing directions to avoid damage.
	Dispose of in accordance with local and national regulations	Follow all local and national regulations when disposing of parts or all of product.
	Temperature Limitation	To avoid damage, store within 14° F to 140° F (-10° C to 60° C).
	Humidity Limitation	Maximum 90% RH (no condensation).
	Temperature Limitation	To avoid damage, operate within 32° F to 104° F (0° C to 40° C).
	Manufacturer	Name of manufacturer.
	Date of Manufacture	Date the product was manufactured.
	Serial Number	Assigned serial number of device(s)
	Catalog Reference Number	Identifies the specific product.

How to properly fit Footbeat

1. The insole pad should be positioned to press into the center of the arch.
2. Adjust the Moc to ensure a firm but not tight fit.


IMPORTANT

A loose or too tight fitting Moc will result in sub-optimal results. Always wear socks while using Footbeat.


How to use and secure the Moc

Footbeat comes with a pair of performance Mocs specifically designed to provide the optimum comfort and performance from Footbeat. The Mocs feature an expanded insole cavity, comfortable toe box and a custom strap that provides closure and support. The adjustable strap enables generous expansion to accommodate swelling. With the strap tightened, the fit should be comfortably firm, not tight, and never painful. The primary purpose of the Mocs is to secure Footbeat in place so that as the insole pad extends into the arch. Footbeat should not push away from the bottom of the foot. The Mocs should be used for minimal indoor walking, while seated or lying down.


1 Charge Insoles

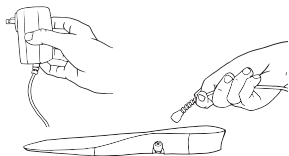2 Insert Insoles into Mocs

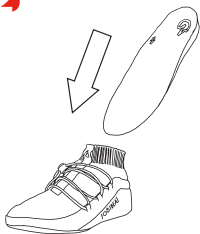
3 Secure Moc to Foot

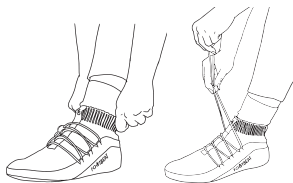4 Turn on with Remote

5 Relax & Recover

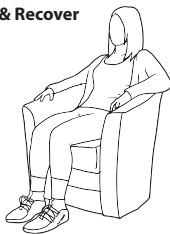Operation


1. INSERT INSOLE into Moc.
2. TURN ON: Push and hold the button on the remote for one second. A green LED light will appear on the remote. Units will activate.
3. TURN OFF: Push and hold the button on the remote for one second. A red LED light will appear on the remote.
4. LOW BATTERY CONDITION: A red flashing LED light will appear on the remote.
5. Footbeat is engaged when you are seated or reclining. Footbeat automatically pauses as you stand or walk. Footbeat resumes as you return to a seated or reclining position.
6. PAIRING: The Footbeat remote works exclusively with the insoles provided to you. However, the remote may un-pair with your insoles from time to time. This can happen due to changing the remote battery or due to a loss of power to the remote. To pair them up again, hold the remote within 6" inches of BOTH insoles, and press the button for one second. The GREEN LED will flash on the remote as with normal TURN ON operation and the insoles will begin to function. If the GREEN LED on the remote continues to flash, move the remote closer to BOTH insoles. Wait up to 1 minute for the Footbeat activate.
7. CHARGING:


IMPORTANT

Only use the supplied charger for charging. For international use, make sure you use the proper outlet converter.


- a. Step 1 – Plug the supplied charger into a compatible electrical outlet.
- b. Step 2 – Align the wand end of the charger to the charger port on the insole. Once the wand is within two inches of the charging port, the magnet will pull the wand into place.
- c. Step 3 – Once the charging indicator light located on the charging port of the insole has turned green the insole is charged. Pull the wand away from the charging port.

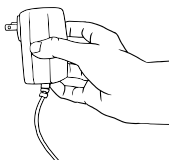STEP 1

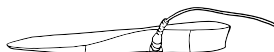STEP 2

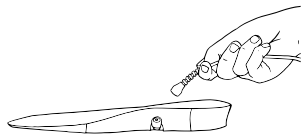STEP 3

IMPORTANT

The rechargeable battery inside the Footbeat is not removable. Attempting to self-service will void the warranty and may cause injury.

NOTE – Achieving a full charge will take approximately 1 hour, which will provide 4+ hours of continual use.

NOTE – Two insoles can be charged at the same time. Due to variations in batteries, the insoles may charge at different rates.

NOTE – Turn off insoles before applying to the charging unit. They can be turned off using the remote.

NOTE – To turn units on or off, the remote must be within 6 feet of the insoles. If more than 6 feet away, the Footbeat may not turn on or off properly.

Caring For Your Footbeat

NOTE – This system does not require any user maintenance. There are no user repairable components. If a problem occurs, please contact Footbeat Support. Opening the case on the control unit, with the exception of replacing the battery, or attempting to self-service any part of the insole will void the warranty and may cause injury.

Routine Cleaning

1. Insole

- a. A damp cloth with mild detergent may be used to wipe the outside of the insoles. DO NOT use a strong cleaning solution that contains bleach. DO NOT put under running water or immerse in water.
- b. To keep the charging port and pad clear of debris, sand and dirt, wipe with a wet cloth.

2. MOC

- a. Machine wash warm. Tumble dry low. Do not bleach. Do not dry clean. Avoid the use of fabric softeners.

IMPORTANT

REMOVE THE INSOLE BEFORE WASHING THE MOC. The insole is not waterproof. Washing the insole will cause the system to malfunction, may cause injury and will void the warranty.

Routine Inspection

Periodically inspect for damage and wear. Check electrical cables and connectors to ensure they are intact and that no conductive wires are exposed. If damage preventing the safe use of Footbeat is evident, please contact Footbeat Support.

Troubleshooting

If your Footbeat stops working, use your remote to turn the power off, then turn the power on to reset.

If reset does not work, then attempt to charge the battery. See charging instructions. Next attempt to turn the power on.

If charging the unit does not work, try replacing your remote battery (CR2032)

If Footbeat is still not working, contact Footbeat Support at 1-844-388-AVEX or email support@footbeat.com.

Frequently asked questions can be viewed at www.footbeat.com

Handling

1. AVOID dropping. Although robustly designed, damage may occur that could cause the system to malfunction.
2. DO NOT open the insole, control unit or charger for any reason.
3. DO NOT submerge the insole in any liquid for any reason. Footbeat has been designed to be water resistant, not waterproof.
4. DO keep the insoles dry.

IMPORTANT

The charging port contains a magnet. Under certain circumstance, magnets may cause interference with some medical devices, including pacemakers and insulin pumps. Keep the charging ports away from such medical devices.

IMPORTANT

The charging port contains a magnet. Under certain circumstance, magnets may cause damage to some electronic devices, including hard drives in laptop computers. Use caution when the charging port is near electronic devices.

Storage

1. Make sure the insoles are turned off and not plugged into a power source.
2. Store in a safe, cool and dry place.
3. This unit contains a Lithium Polymer Battery. Please observe the following Battery Storage Cautions:
 - a. Do not store packs in places of high temperature or under direct sunlight.
 - b. Do not store in place which may expose them to rain, water or high humidity.
4. The Footbeat system (insoles, remote, socks, Mocs, Quick Start Guide) should be stored in the provided tote bag to minimize potential damage.
5. DO NOT store the system in a car where the temperature could exceed the recommend storage conditions. This could damage the units and void the warranty.

Refer to the technical information and specifications section for proper storage conditions.

Disposal

1. When the device has reached the end of its useful life, please dispose of properly in accordance with all local and national regulations.

Standard Features	Auto-adjusts to arch height Auto-OFF when walking Insoles contain rechargeable lithium ion battery Estimated 4 hours of use between charges Cordless device when in use Remote communication via Bluetooth (low energy) Remote on/off with CR2032 replaceable battery
Mechanical Dimensions Controller Charger Moc with Footbeat (each) Weight	1.25"W x .25"D x 4.5"H (3.2cm W x 0.6cm D x 11.4cm H) 1.75"W x 3"D x 1.5"H (4.4cm W x 7.4cm D x 3.5cm H) 4"W x 6"H (10.2cm W x 15.2cm H); Length varies on size 1.61 lb (0.73 kg) net, while in use 4 lb (1.81 kg) net, packaged box
Environment Operating Storage	32° F to 104° F (0° C to 40° C); Maximum 90% RH (no condensation) 14° F to 140° F (-10° C to 60° C); Maximum 90% RH (no condensation)
Power Input Voltage Input Frequency Input Current Output Voltage Output Current Power Consumption	(NOTE: Values below pertain to recharging cycle) 90-264 VAC 47-63 Hz 0.6 A 9 VDC 0.66 A 6.0 W

Refer to EMC Declaration included with your Footbeat

FCC Part 15 Information to User

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and

(2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by AVEX could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

*Reorient or relocate the receiving antenna.

*Increase the separation between the equipment and receiver.

*Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

*Consult the dealer or an experienced radio/TV technician for help.

This portable transmitter with its antenna complies with FCC/IC RF exposure limits for general population / uncontrolled exposure.

Limited Warranty

Footbeat must be under warranty to qualify for replacement.

30-Day Full Replacement Guarantee

If you are not completely satisfied with Footbeat you may return the complete system (two insoles, remote, one pair of Footbeat Mocs, socks, charger and tote bag) in the original packaging, no questions asked within 30 days for a replacement set.

EXTENT OF LIMITED WARRANTY - This limited warranty of the Footbeat does not cover damages due to external causes, including without limitation misuse, alteration or repair by a non-authorized person; accident or use not in accordance with the product instructions as well as failure to follow the care and cleaning instructions in the manual. AVEX warranty does not cover damage due to general wear and tear, burn marks, punctures, cuts, and abrasions from outside causes or misuse and/or abuse by the wearer.

To obtain warranty service please do one of the following: go to www.footbeat.com/support and following the instructions, contact Footbeat Support at 1-844-388-AVEX, or email our Footbeat Support at support@footbeat.com.

THESE WARRANTIES REPLACE ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. AVEX MAKES NO EXPRESS WARRANTIES BEYOND THOSE STATED HERE. AVEX DISCLAIMS ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. IF YOU RESIDE IN A JURISDICTION THAT DOES NOT ALLOW THE EXCLUSION OF IMPLIED WARRANTIES THEN THIS LIMITATION MAY NOT APPLY TO YOU.

THIS WARRANTY IS ONLY EXTENDED TO THE ORIGINAL OWNER AND IS NOT TRANSFERABLE. ORIGINAL PROOF OF PURCHASE IS REQUIRED BEFORE A WARRANTY CLAIM MAY BE PROCESSED. WHEN AVEX PRODUCTS ARE PURCHASED FROM UNAUTHORIZED RESELLERS OR DISTRIBUTORS, THE PURCHASED PRODUCT SHOULD BE CONSIDERED AS IS AND WILL VOID THE WARRANTY.

Footbeat Support

Your satisfaction with Footbeat is important to us. Please feel free to call and talk to a Footbeat Support staff member with questions or comments. We cannot offer an opinion or propose a recommendation about your individual health situation. Please consult your healthcare provider, athletic trainer, team coach or physiologist for specific guidance about your situation.

AVEX LLC., Grand Junction, CO • www.footbeat.com

Footbeat Support

Toll Free: 1-844-388-AVEX

Hours: 8:00 A.M. – 5:00 P.M. MST

Email: support@footbeat.com

Document Number : 600-0001-01 (Rev. A)

© Copyright 2017 Avex, LLC. All rights reserved.

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Documento número: 600-0001-01 (Rev. A) © Copyright 2017 Avex, LLC.